

simple&smart

ANYA 

Cuando menos es más



SIMPLE inside
Less assistance need
easy installation
LESS electronics
for a **MORE RELIABLE UNIT**
That's why we are
simple&smart





1 CHARACTERISTICS



CHARACTERISTICS

COMFY & COZY

improving ergonomics

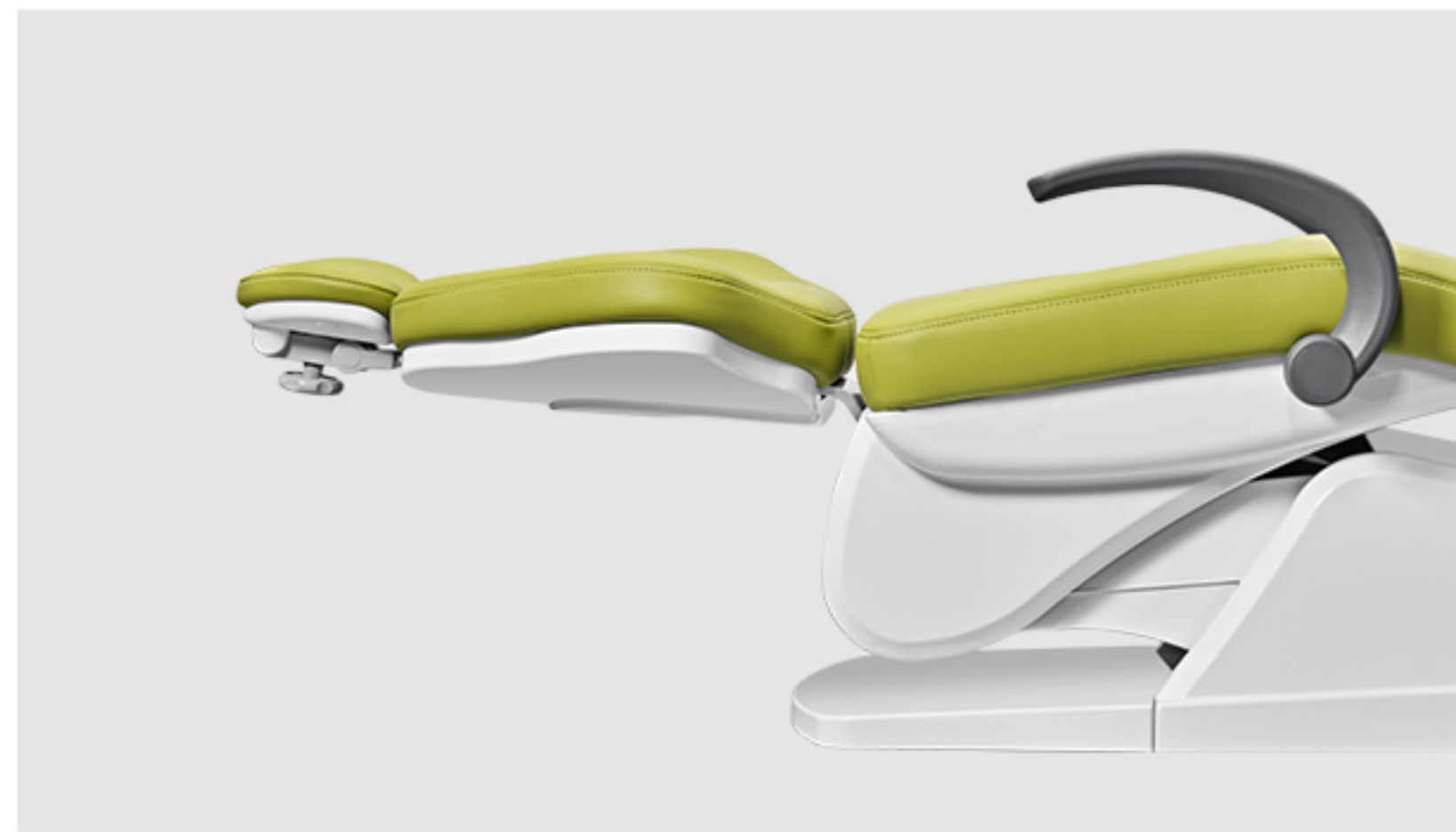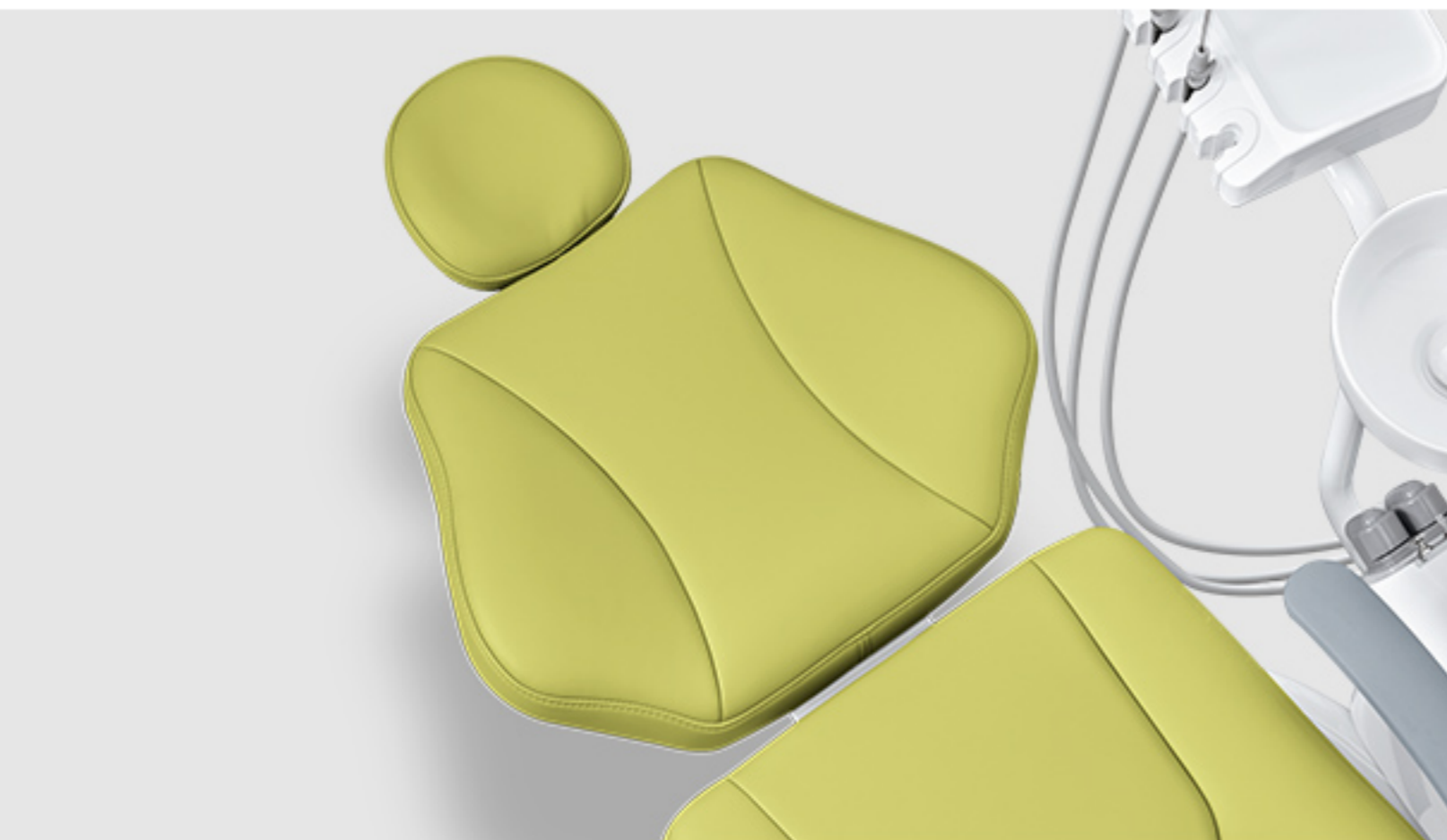

The **SEAT** and the **BACKREST** have been designed to improve the patient's comfort and stability on the chair. The particular ergonomics allows the doctor to work in a more suitable position. Working in 12 o'clock position will be easier

The structure of the **BACKREST** is extra Sturdy design improves the ergonomics of the seat and makes the approach of the


doctor to the patient much more easier, in this way the doctor can reach a **better working position**.

The design of the **CHAIR BASE** grants a wider movement of the seat. It is thin and allows the chair to reach a lower position, to grant a wider working space and improve the doctor 's comfort during the daily practice.

**CHAR
ACTER
ISTICS**

GOOD TO THE BONE

resistent structure, performing features

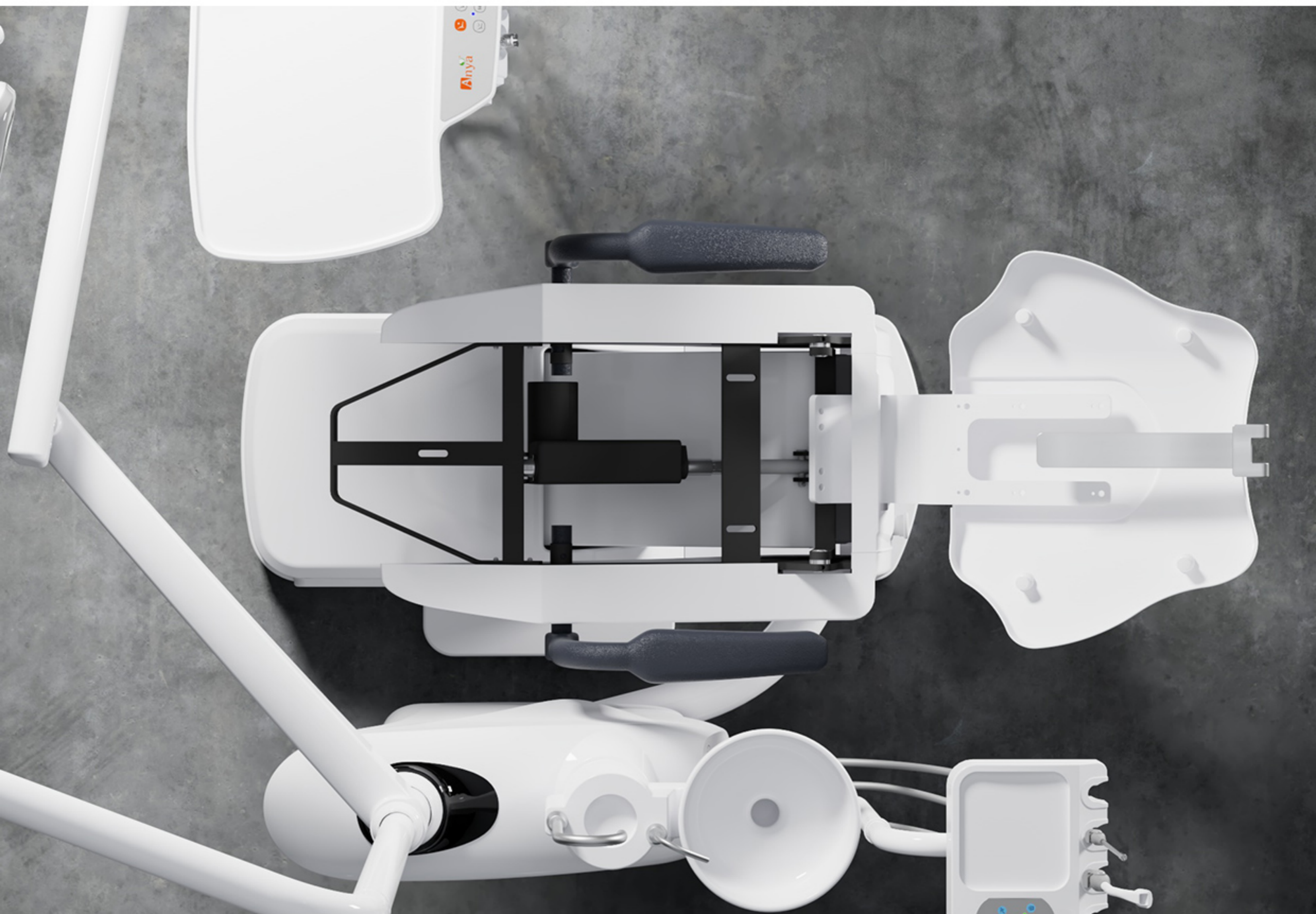

The Dental Unit A1000 Is Equipped With Two Separated Airtight Filters,one For Each Suction Cannula, Which Grant A Perfect Suction Action And Hygie-nic Standards.New, More Powerful And Less Noisy Motor

CHARACTERISTICS

WHAT'S DIFFERENT

guess what's new

The A1000 has been designed to meet the needs of those customers who need very **SIMPLE** and **FUNCTIONAL** devices, essential but at the same time **EFFICIENT** and **LONG-LASTING**. Designed to further contain costs, A1000 is **PERFECT FOR TENDERS** where excellent quality standards are needed besides excellent options and advantageous prices. So, what differentiates A1000eco from other brand units?



POWER SWITCH

Unlike the other brand units, the A1000 does not have the lever to activate the air (from the compressor) at the chair base, it simply has a switch for turning ON and OFF the power supply. Also easy to switch to Junction Box.

TRAY ARMS

Both the instrument tray and assistant tray arms have been modified to simplify the structure of the unit. A **smoother desing** and a **simplicity in movements** are the main features of this unit. The lamp pole is mounted over the instrument tray arm.




CHARACTERISTICS

WHAT'S DIFFERENT

guess what's new

Rotatable Cabinet

Can rotate outwards, convince to dentist and assistant, coordinate with each other during treatment.

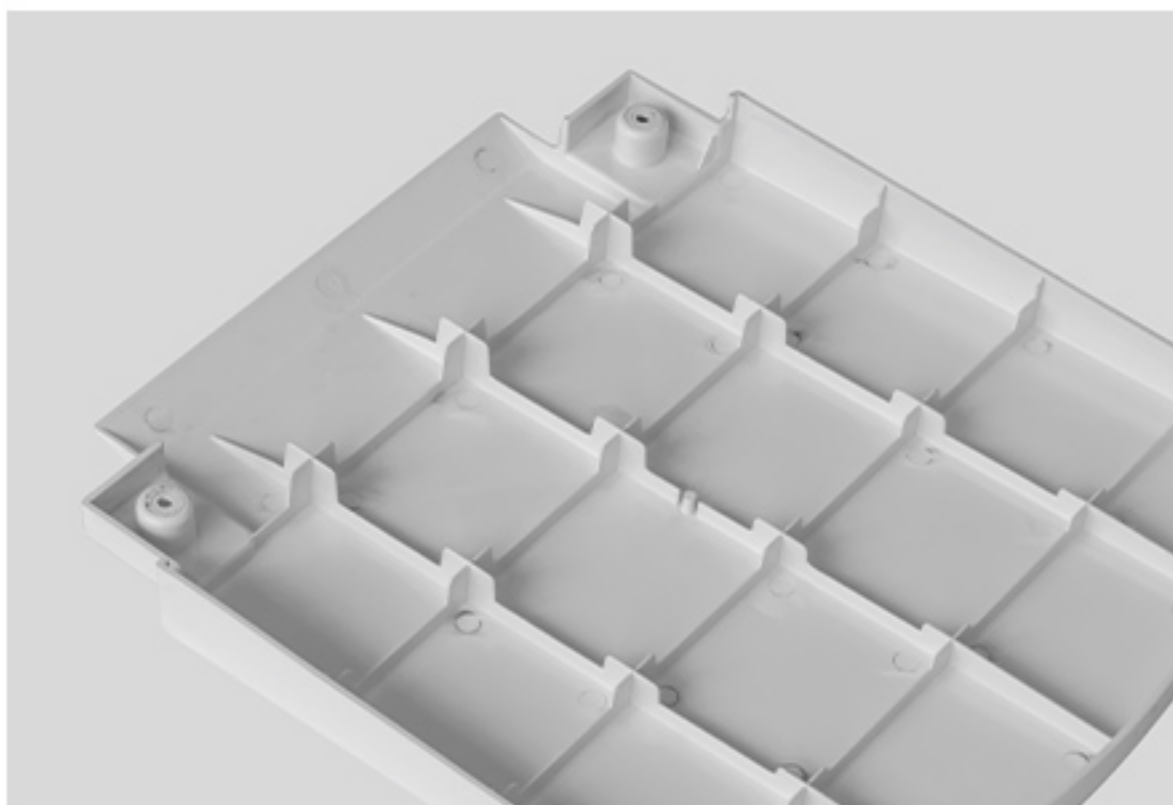

All ASA Plastic By Injection Moulding

Compared with ABS, Long life-time, Acid and alkali resistance, Cost double time than ABS Material

Self-Designed PCB System Import Chip

Own-Designed control system which problems can be checked easily, transmission stability and fast so system not easy to crash



CONF IGUR ATION



CONF
IGUR
ATION

FIT YOUR SPACE

pick what you need

You can choose which kind of **INSTRUMENT TRAY** is more suitable for you, based on your habits. Also the choice of the **SPITTOON BOWL** is up to you: the dark glass spittoon bowl is removable to simplify the cleaning operation.

HANGING TYPE



TOP MOUNTED




CERAMIC SPITTOON



GLASS SPITTOON

CHARACTERISTICS

STANDARD

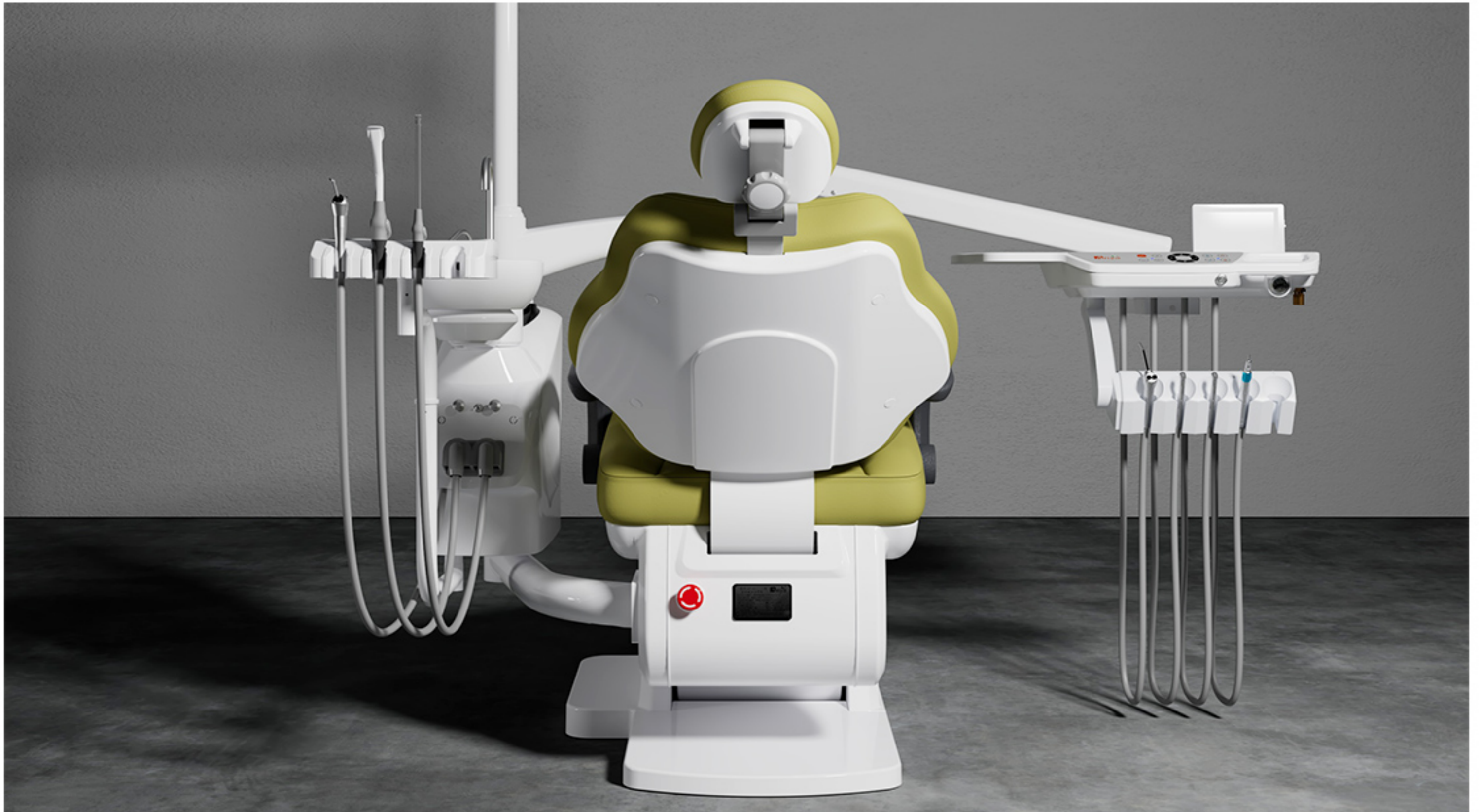

Standard Equipment

- Strong & Weak Suctions·1 pc each
- Multifunctional pedal switch·1 set
- High speed handpiece tube·2 set
- Low speed handpiece tube· 1 set
- High class dentist stool· 1 set
- Dental chair auto spittoon flush system·1 set
- Cup filler control constant temperature system·1 set
- Assistant control system·1 set
- Three-way syringes·1 pc
- LED flm viewer·1 set
- Luxury LED reflective sensor light·1 set
- Purified water supply system·1 set
- 90°ceramic rotatable spittoon·1 set

CUS
TOM
I Z E

your
dental
unit



CUS
TOM
I Z E

EXPRESS YOURSELF

you ask - we do



Our dental units can be equipped, according to your needs, with the **addition of many options**.

The internal simplicity of the machine remains the same but, following your requests, **we can customize your unit in order to make it perfect for you!**

Your experience and feedbacks help us to constantly improve our unit and find the best solution in any situation

We are the factory which can fulfill these requests that lie outside the standard production process. Therefore it is possible to adjust ergonomic and functional aspect of the unit according to your requests

If you have a wish, we try to make it real!

OPT ION ALS

accessories

- lamps
- media kit
- handpieces
- micromotors
- valves
- instruments cart
- maintenance

FIBER LEATHER

More comfort with seams! The standard cushions in PU Leather are manufactured without seams as they can be cleaned more easily. On the other hand, the Fiber Leather seat grants an **higher sense of comfort** to your patients.



EXTRA STOOL

An extra stool can be suitable for your assistant. Different colors are available to meet your requests and your needs

OPT ION ALS

accessories
lamps
media kit
handpieces
micromotors
valves
instruments cart
maintenance

ECO LED PEBBLE

Equipped with sensor for intensity regulation. The intensity of this lamp can reach 40.000lux. Easy to clean and with a pleasant design, this lamp is the perfect example of excellent value for money.



LED LAMP EVA

High quality LED lamp - Illuminating power from 3.000 to 50.000 lux. You can select the ideal color temperature for dental practice 5700k, 5000k, 4000k, and 2700k. The patented, unique, reflected-light optical design ensures a well-defined light spot, preventing glare for patients. Its 3D head rotation >360 allows to illuminate every corner of the mouth.



LED LAMP MAIA

High quality LED lamp, from 8.000 to 45000(lux. Maia is the ideal choice for all dentists searching for a professional LED lamp, a reliable solution, long lasting and with high-level lighting performance at a competitive price. The accurate definition of the MAIA dental lamp luminous spot eliminates the risk of glare, safeguarding the eyes of patients and medical team.



OPTIONALS

accessories
lamps
media kit
handpieces
micromotors
valves
instruments cart
maintenance

ENDOSCOPE KIT

- **Monitor** 24 V, color, 15". with glass protection.
- **Intraoral USB Camera:** Full HD resolution, auto focus and mouse, LED illumination -by simple&smart
- **Arms Set** for monitor installation


MONITOR and CAMERA are installed on the dental unit with a special pole set .The camera allows to save images on a PC or into the dock-station supplied with the device,thanks to its USB connection cable. This Full HD camera (1920 x1080) is equipped with autofocus and mouse,LED illumination with 3 different intensity levels.The head of the camera can rotate 280°.this feature facilitates your job and reduces discomfort for the patient.The Endoscope KIT is the PERFECT SOLUTION to show to your patients your work in real time.



SUPERMEDIA KIT



- **Panel PC** Installed on the unit, touch screen, wire-less mouse and keyboard.
- **USB Camera:** Connected to the panel PC
- **Arms Set** for monitor installation

MONITOR and CAMERA can be installed also se-parately, but the BEST SOLUTION is to use them together. The PC is touch screen and it works as a normal computer. Internet connection,RVG system,endocamera and all management programs can be controlled directly from your operation site. The combination with the camera allows to get and save all images you need.

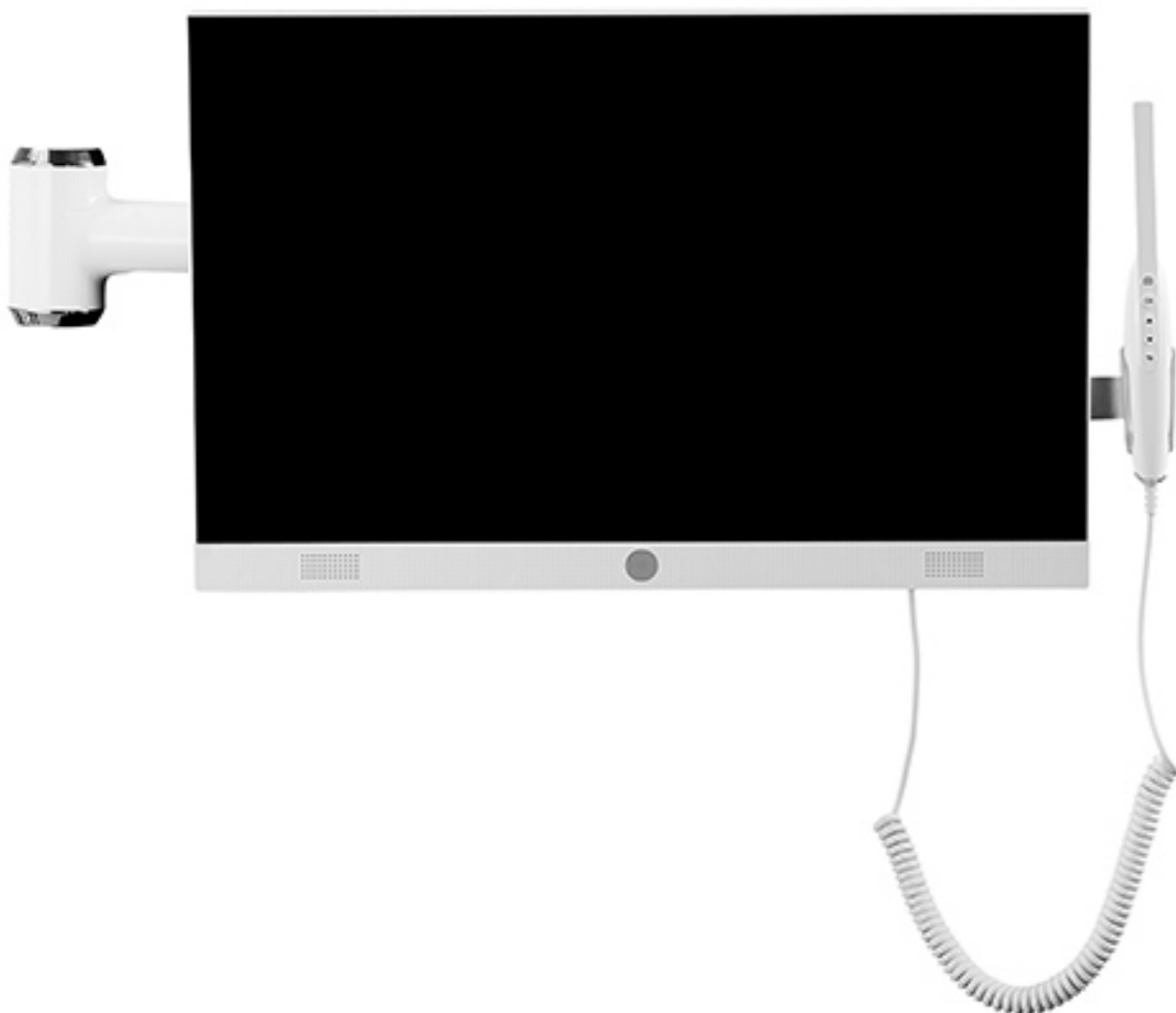

OPT ION ALS

accessories
lamps
media kit
handpieces
micromotors
valves
instruments cart
maintenance

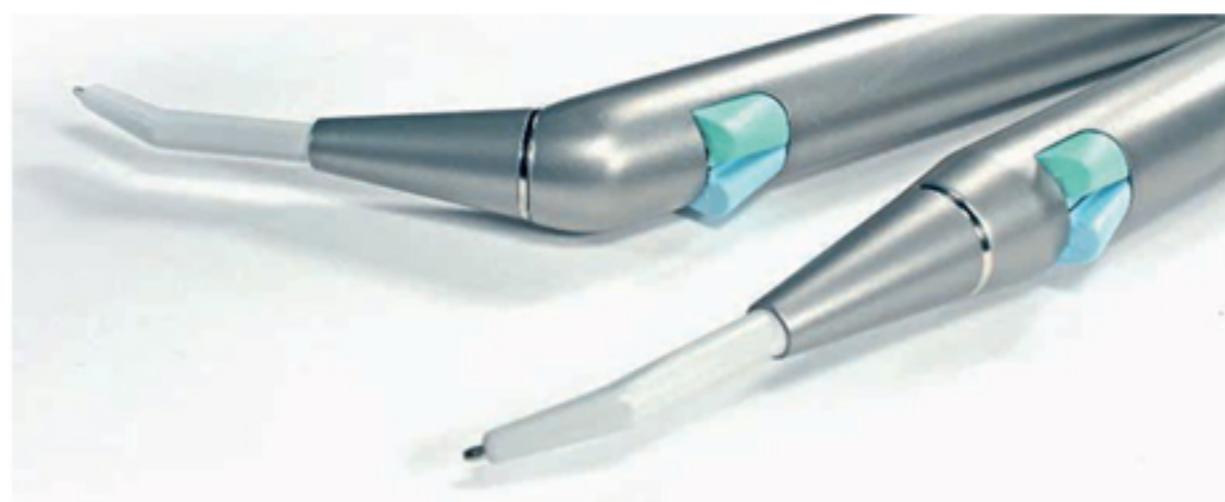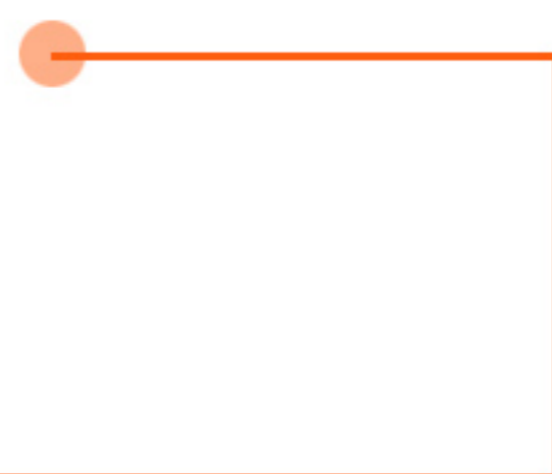

ENDOSCOPE KIT

Straight or curved
3 or 6 functions



High power



Safe and clean



Easy to grip



Low noise




Solid quality



Convenient to clean

OPT ION ALS

accessories
lamps
media kit
handpieces
micromotors
valves
instruments cart
maintenance

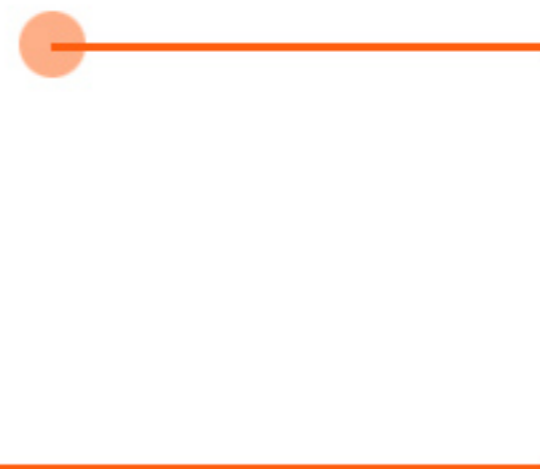

WOODPECKER **SCALER**

with light
EMS compatible
Satelec compatible on request



WOODPECKER **SCALER**

without light
EMS compatible
Satelec compatible on request



WOODPECKER **CURING LIGHT**

With cable on assistant tray



ANYA **CURING LIGHT**

Second Mode, Normal Mode,
Caries Mode, ulse Mode, ulse Mode,
Intelligent Main Switch

OPT ION ALS

accessories
lamps
media kit
handpieces
micromotors
valves
instruments cart
maintenance

On all our dental units we install **ANYA** valves, according to the specific requests of each customer and to the type of suction system present in the clinic. Our dental units are preset by default to be connected to a liquid or wet suction machine. In addition, it is possible to integrate the unit with.




1. PLACE SELECTION VALVE

In case of multiple installations (two units or more connected to a single suction motor).

2. BUILT-IN AMALGAM SEPARATOR

Built-in amalgam separator, elegant and efficient 97.9% separation rate, meet EU standards & requirements. Convenient for operation and maintenance to lock, simple and easy to operate.

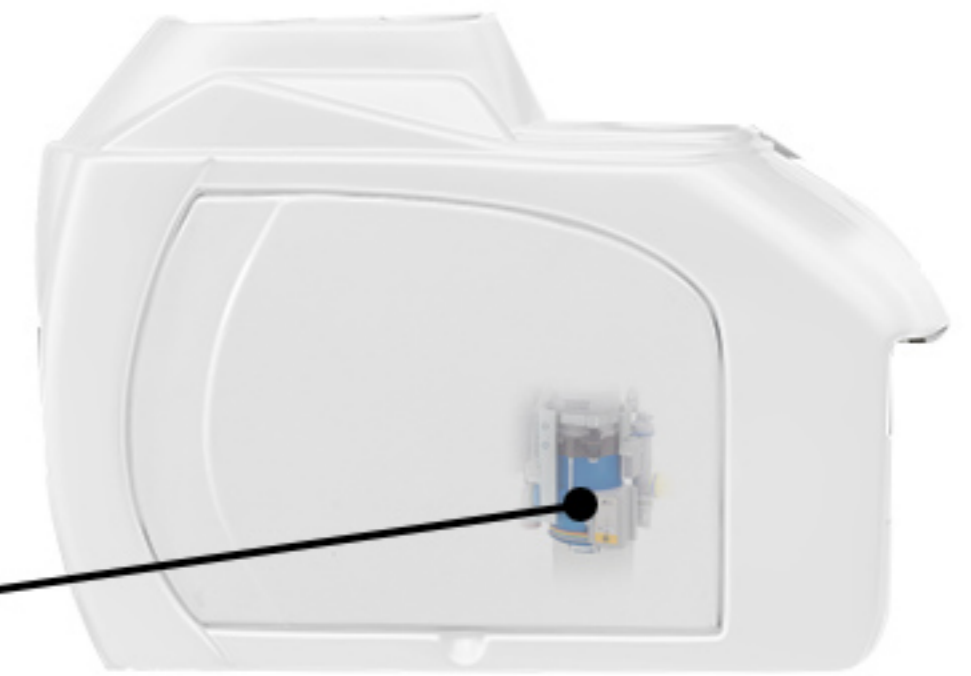

3. SPITTOON COLLECT VALVE

Conveniently collect dirt from the spittoon and better maintain the dental chair

It is very important to analyze the situation of the clinic before ordering the unit. **The ANYA staff will help you and suggest you the best solution to fit the unit to your system.**

OPT ION ALS

accessories
lamps
media kit
handpieces
micromotors
valves

instruments cart
maintenance



In replacement to the instrument tray you can have your dental unit with the instrument cart The cart has the same equipment as the instrument tray:

• FIE CORDON MODULES ON DOCTOR TRAY:

- 1 syringe with 3 functions
- 1 woodpecker scaler, EMS compatible with 5 tips included
- 3 cordons (2 with optic fibre, 1 without optic fibre) midwest type

• 15 FUNCTION CONTROL PANEL

• X-RAY VIEWER


OPTIONALS

accessories
lamps
media kit
handpieces
micromotors
valves
instruments cart
maintenance

Our aim is to be as fast **as possible** **in assisting you**. That's why we have all these spare parts ready for the delivery and they can be shipped within 12 hours from the reception of the request.


In order to be even faster in solving your possible problems, we created **2 SPARE PARTS KITS**. If you keep these kits in your warehouse or clinic, you can solve any possible problem in a short time.

BASIC SPARE PARTS KIT

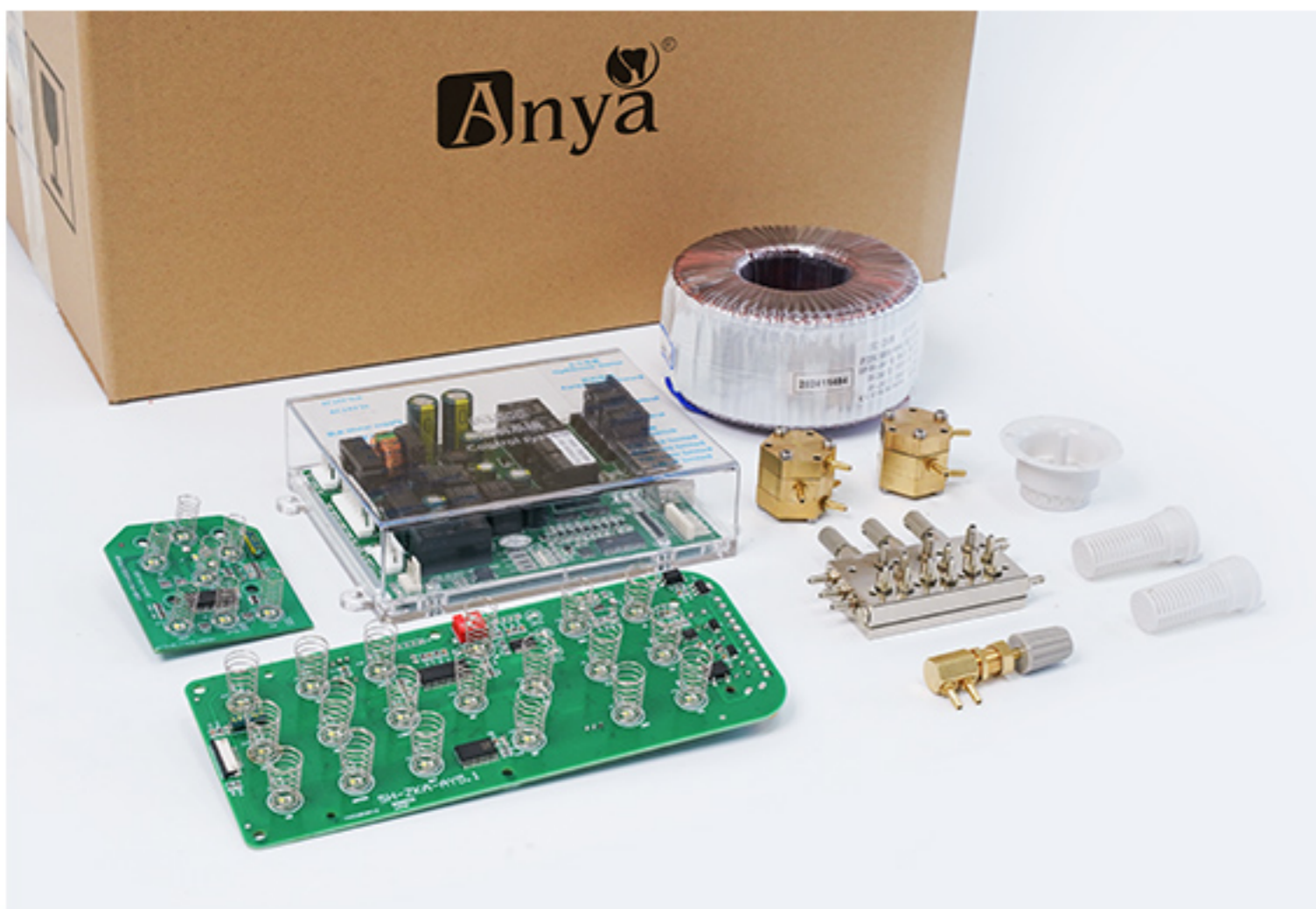


- All the electrical boards of the unit (BASIC SPARE PARTS KIT)
- 1 Main Board
- 1 Transformer
- 1 Assistant Tray PCB
- 1 Instrument tray PCB
- 2 Water Valve
- 1 Triplex Diaphragm Valve
- 2 Suction Filters
- 1 Basin Filter
- 1 Regulation Valve

COMPLETE SPARE PARTS KIT

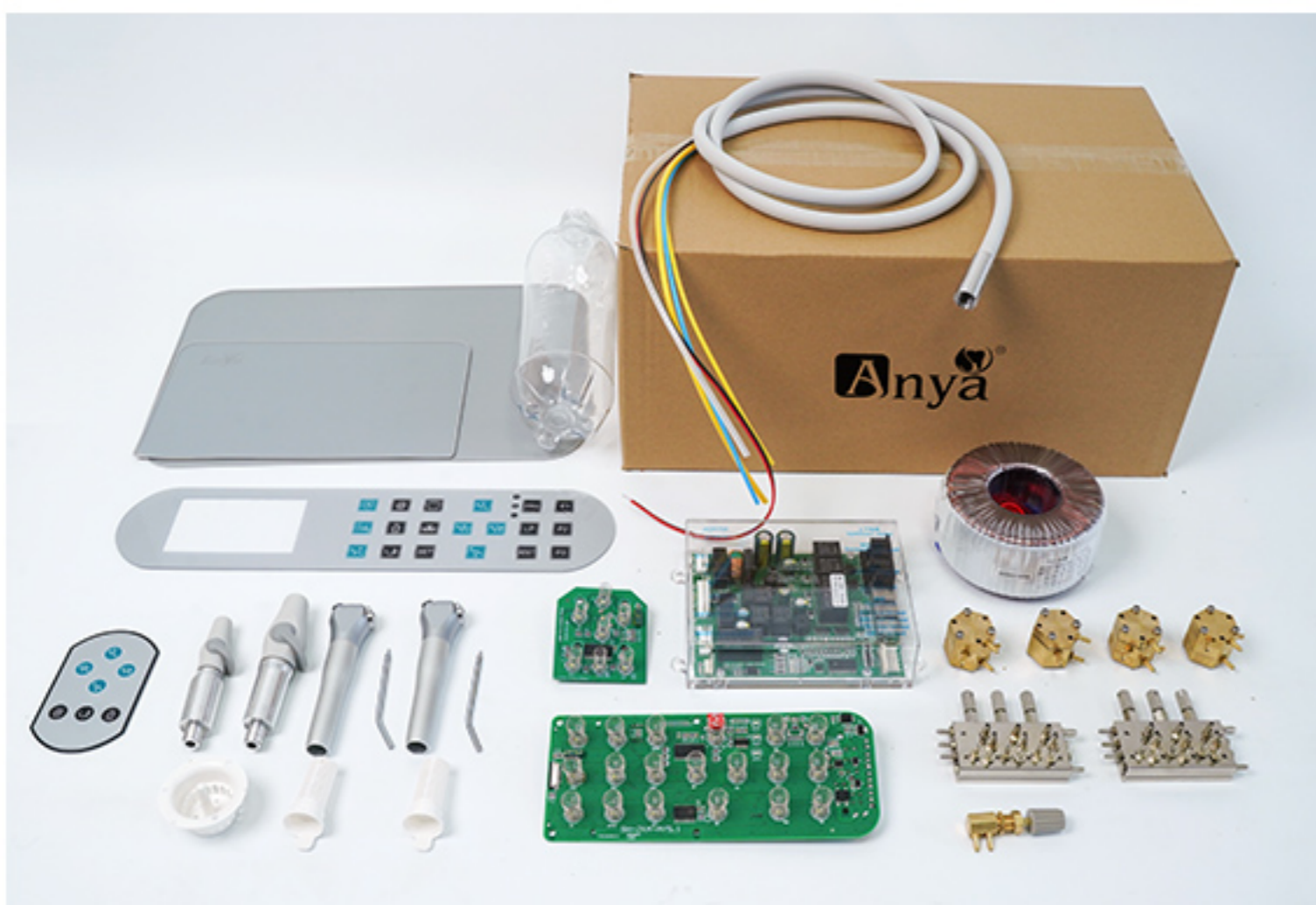

- 1 Mainboard PCB
- 1 Transformer
- 1 Assistant Tray PCB
- 1 Instrument Tray PCB
- 4 Water Valve
- 2 Triplex Diaphragm Valve
- 2 Suction Filters
- 1 Basin Filter
- 1 Regulation Valve
- 1 Weak Suction Handle
- 1 Strong Suction Handle
- 2 Syringe
- 1 Optical Fiber Tube
- 1 Assistant Tray Label
- 1 Instrument Tray Label
- 1 Bottle For Disinfection System
- 1 Set Of Silicone Pads
- 1 Tap For Basin
- 1 Tap For Glass

4 COLOR YOUR CLINIC

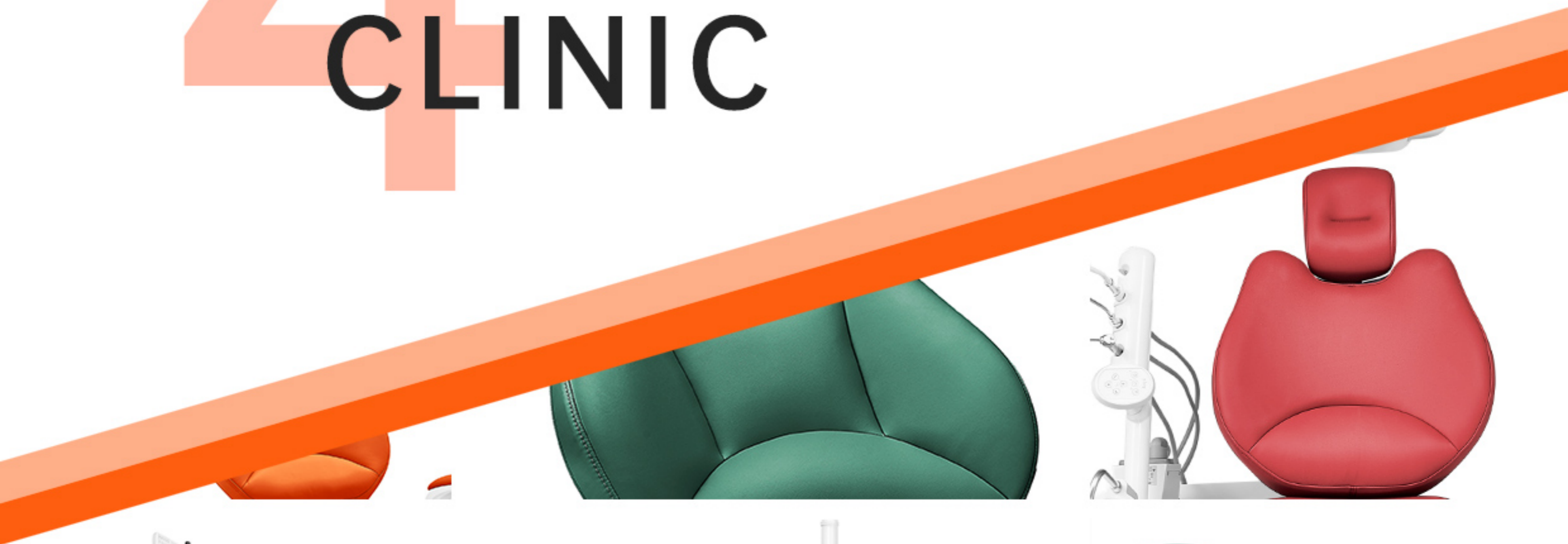

COLOR YOUR CLINIC


IT'S A MATTER OF LEATHER

seams or no seams, that is the question!

Standard PU, standard Fiber Leather and more than 200 customized colors: ANYA has the widest choice of upholstery you have ever seen!

The Pu standard polyurethane upholstery is available with no additional costs. The standard PU leather is seamless for an easier cleaning procedure. With a small additional cost you can select the elegant Fiber Leather upholstery, with seams to grant to your patients a higher comfort.

PU LEATHER

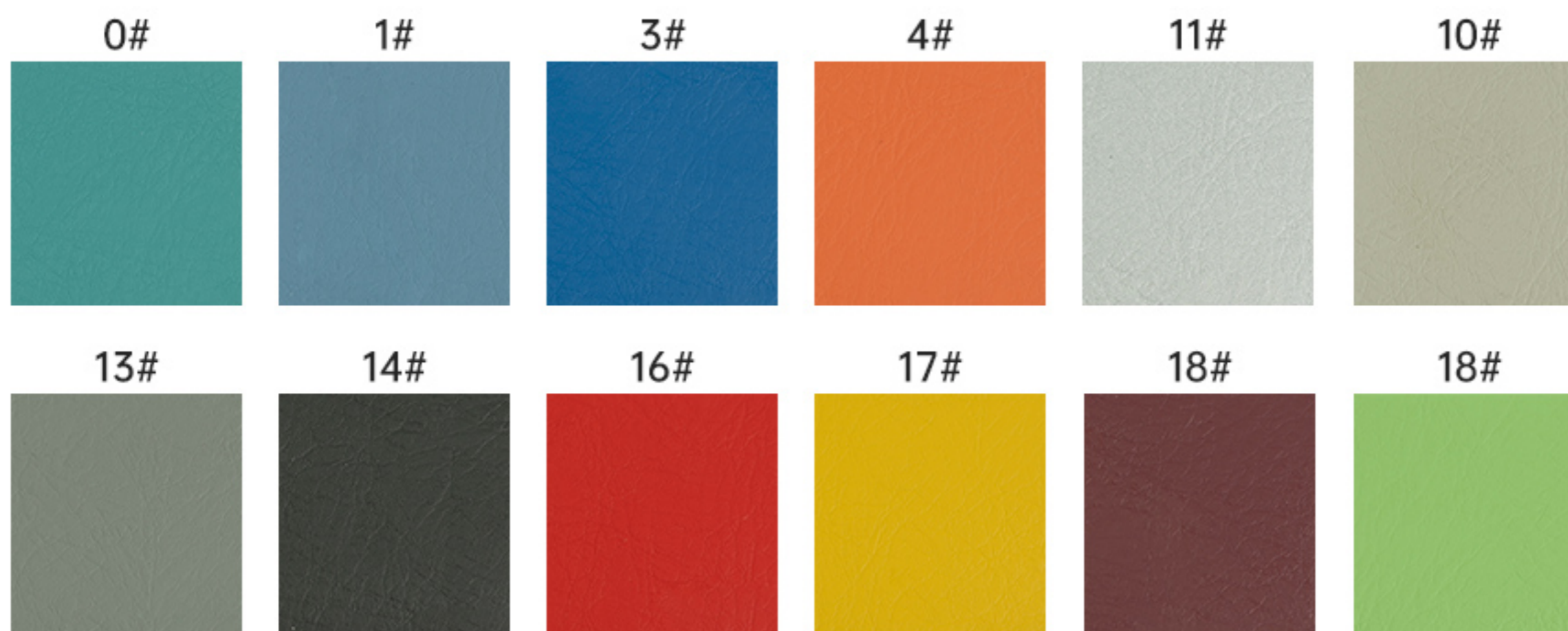

FIBER LEATHER



5 COMPANY



COM
PA
NY

012-72 AY-A48001
ASR2412359

012-72 AY-A48007
ASR2412368

012-64 AY-A48001
ASR2412246

007-11 AY-A48007
ASR2412248

禁力叉车



AT YOUR SERVICE

PRODUCTION

The production process is the heart of our company, and for this reason we follow every stage with particular care, starting from the choice of materials, paints and components, to assembly and final control. Each unit, in fact, is subjected to a careful double check, carried out by our specialized technicians, during assembly phase and before shipment. Our ISO 13485 standards testify that these processes are actually carried out and respected.



STOCK

We want to be extremely fast in deliveries, so we always have in stock several units, in large part already assembled and ready to be configured according to our customers' needs. All the necessary parts and all the options are available in stock and we can complete the assembly of the unit and arrange the delivery in a few days.

SHIPMENT

The shipment can be arranged according to each customer's needs. We usually ship the unit in a single wooden case, but, if required, we can ship the unit in two separated boxes in order to make the delivery to the clinic easier. The two separated wooden cases are smaller in size and less heavy, so the download from the truck is easy even without a lifter and the reduced dimensions allow an easy storage while waiting for the installation.





ASSISTANCE

Always at your side

WARRANTY EXTENSION

The standard warranty on our dental unit covers a period of 1 year but, if you want you can extend the warranty up to 5 years

WATCH OUR VIDEOS

On our website and on our YouTube page you can find many videos of the unit, on its configuration, on the options available, tutorial videos for the installation, maintenance, sanitization, etc., in order to set up and take care of your chair in the best way possible

WWW.anyadent.com

We created this website to support our dealers and customers in the best way possible. You can have a complete overview on our company and our products you can even configure your customized unit. Using our configurator you can receive a personalized quotation without any purchase obligation.

With Simple&Smart you can always find the best and simplest solution for your needs.

Follow our social media pages to discover more and to be updated on our news, products and events.

LIVE ASSISTANCE

WhatsApp number on request if you are our dealer

REPORTS ON FILES

Complete PDF reports are available to be sent by email. Detailed information for problem solving are explained in a specific manner on the bases of the problem you have indicated





CERTIFICATES

we take it seriously

ISO & CE CERTIFICATES

Our dental units are Class IIa medical devices. The **EC -07002 Certificate** and **ISO 13485 Certificate** ensure that our products and procedures are compliant to the regulations in force and to these severe ISO standards. Our dental units are certified CE1370. These certifications guarantee the quality of our products and services

 **Anya** Foshan Anya Medical Technology Co.,Ltd

EXPORTDEPT:
+86-757-81803709

WEBSITE:
www.anyadent.com

E-MAIL:
adam@anyadent.com

ADD:
Floor 3 and 4, Building 9, South of Hengyi Road, North of Hengshisi Trunk Road, Nanzhuang Town, Chancheng District, Foshan City, Guangdong Province, China

